

How To Use Your Oticon Intent hearing aids at the Fox Theatre in Atlanta


The Fox Theatre in Atlanta has an Auracast™ audio system that allows you to walk into the venue and connect instantly to stream the performance audio directly through your hearing aids for crystal-clear sound.

Requirements:


1. Oticon Intent hearing aids
2. Your hearing aids must be paired with your smartphone using the Oticon Companion app.

Instructions:

1. Open your Oticon Companion app
2. Tap on the ●●● three horizontal dots on the bottom right

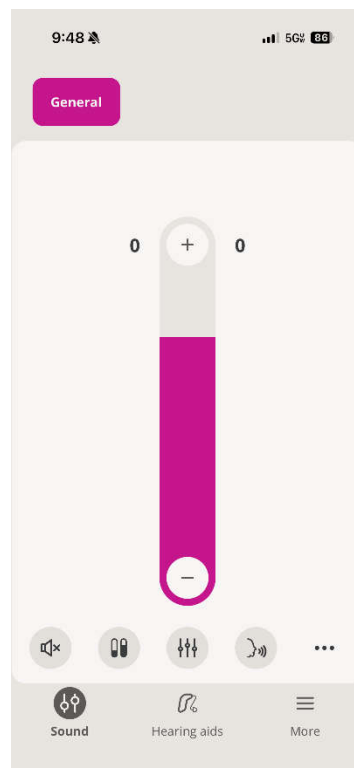


3. Select the **Auracast Broadcast** option at the bottom

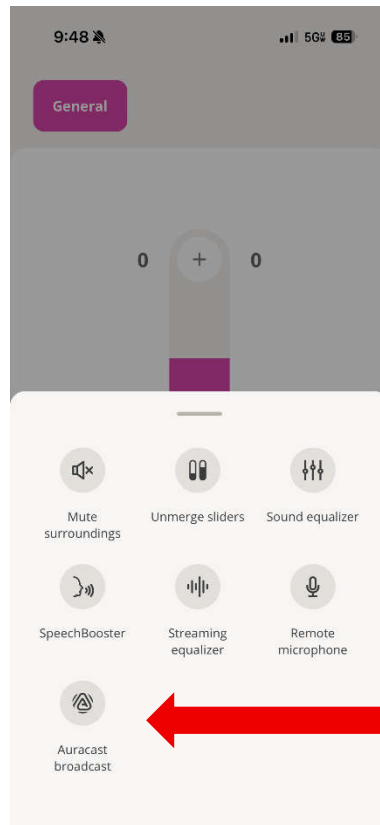


4. Search for broadcasts nearby

*Note: Some venues will have a QR code to access a broadcast

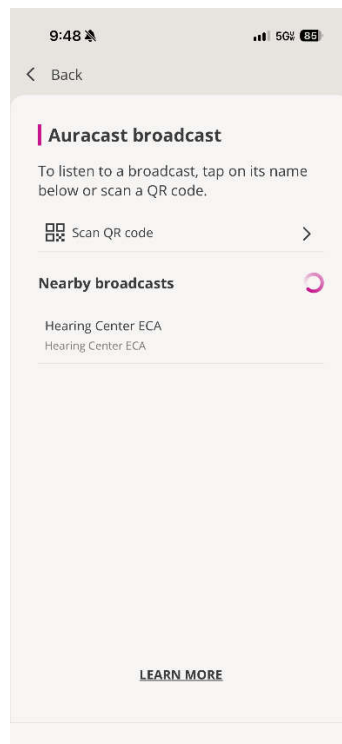

5. You can **adjust the volume of the Auracast broadcast** on the volume control.

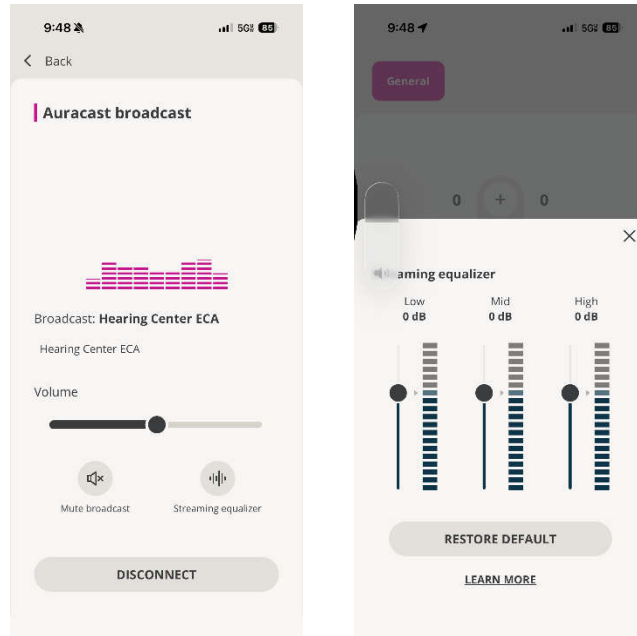

You can **adjust the volume of your surroundings** which changes the volume of your hearing aid microphones on top of your ears on your normal volume control on the home screen (ex. someone behind you is eating popcorn or talking and you want to turn them down).

