



HEARING CENTER


Eye Consultants of Atlanta

How to disconnect then re-connect your hearing aids from your iPhone at home


Step 1: Go to Settings App




Step 2: Go to Accessibility

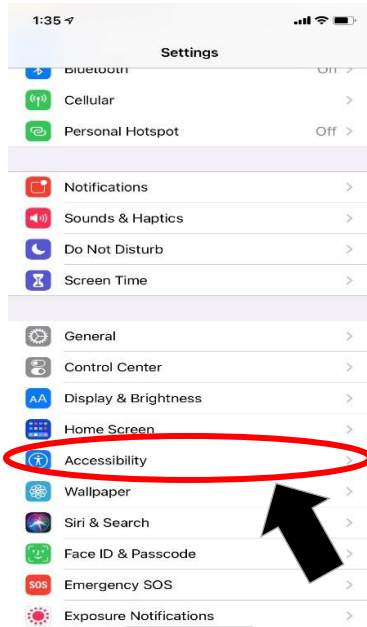


Step 3: Go to Hearing Devices

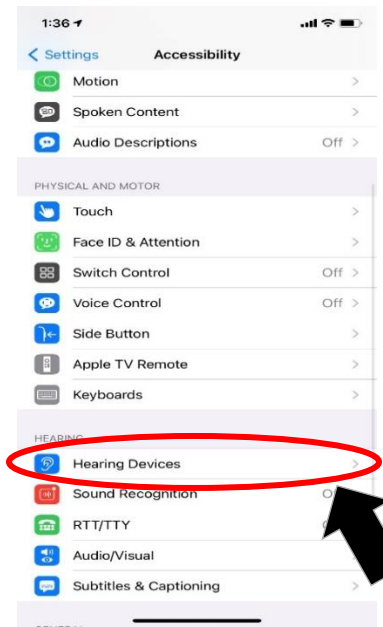


Step 4: Tap on your name

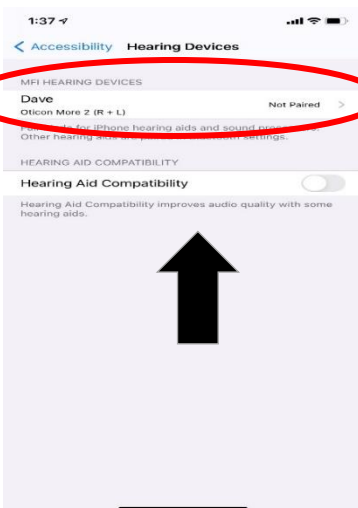

Step 5: Scroll down and hit "Forget this device", scroll back up then hit the "< Hearing devices" in top left to go back a screen

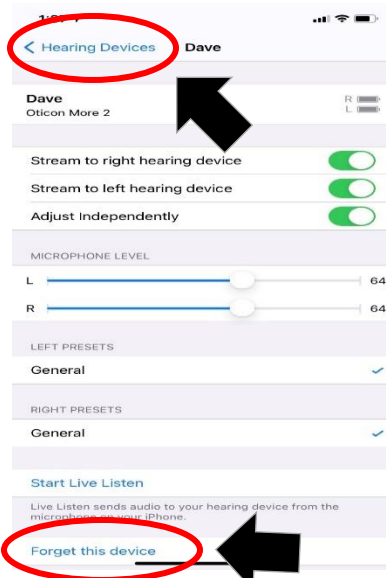


Step 6: Turn off your hearing aids by pressing buttons on BOTH aids SIMULTANEOUSLY for 5 seconds, then turn BOTH back on by pressing BOTH buttons another 5 seconds. If aids are non-rechargeable, open & close battery doors instead of pushing buttons.

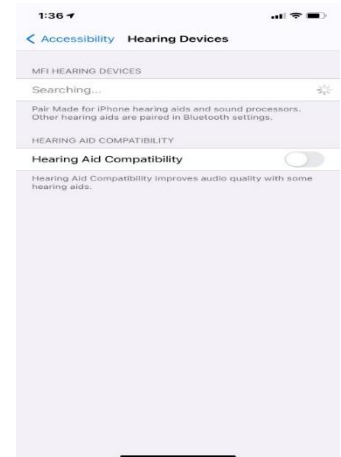


Turn page over




Step 7: Your name should appear. Once you see (R+L), then tap on your name.

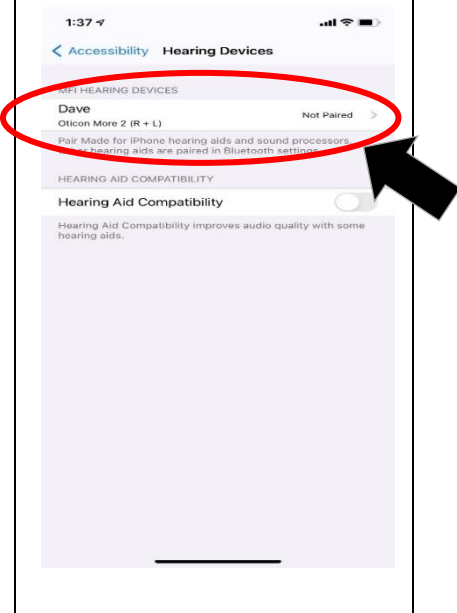


Step 8: Tap on Pair, once for each hearing aid

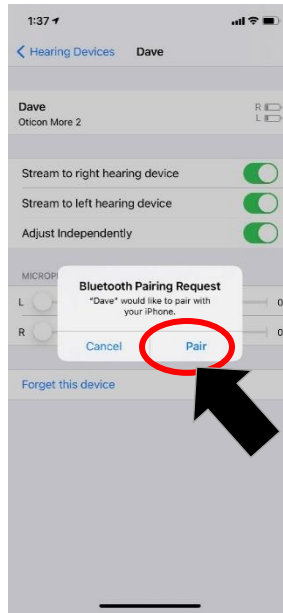

Step 9: Push "< hearing devices", now that aids are connected

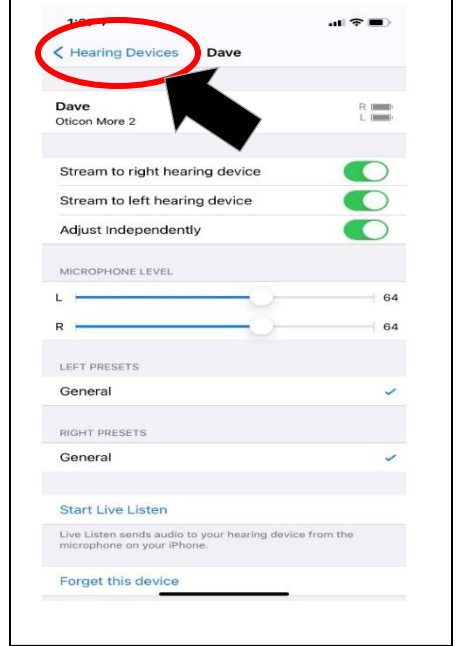

Step 10: You will see it says both aids are connected

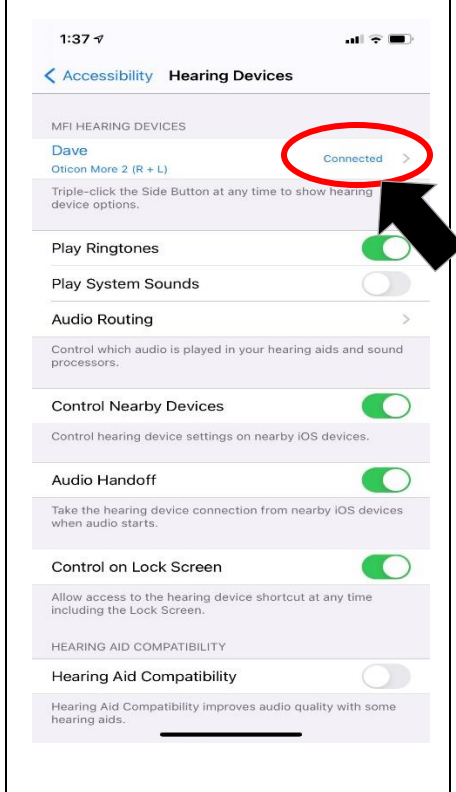

TIP: Turn off play ringtones so that when the phone rings, it rings from your phone and NOT in your hearing aids. This allows you to be able to find your phone and is less likely to frighten you!

TIP: Leave Hearing aid compatibility OFF!

